


Enlarge word tile with click

Step-by-step instructions for Enlarge word tile with click in Express Me.

Watch the video, follow the steps, or open the complete transcript for a text version of the instructions.

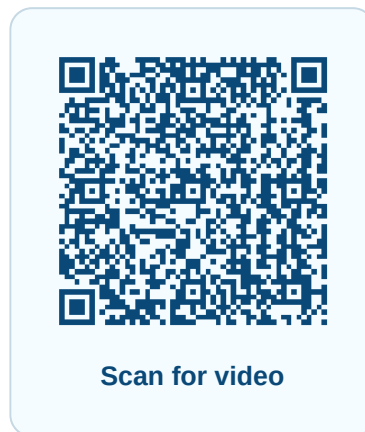

Step-by-step instructions

1

Open the app settings

Tap the blue profile icon in the top right and choose **App settings**.

2

Turn on enlargement

Scroll to **Enlarge image on touch** and set the switch to **On**.

3

Test it

Go back and tap a word tile; the image is now enlarged.

Complete transcript

▼ Read the complete video transcript

Hi, I'm going to show how you can set a word tile to enlarge at the moment that you click on it.

Right now, when I click on it, it's only spoken. Yes.

But we can also go to the settings by clicking on the blue figure at the top right, then to app settings. If we scroll down a little bit, you see here: enlarge image when touched.

When we turn that on and then go back by clicking the X at the top right, we can now tap on yes. Yes.

And then you see that the tile is enlarged. This makes it easier to actually physically show the word you want to show to the person you're communicating with.